

User Manual


Mandala Wand GEN.2



Please read the manual before using the tools.

Version 1.0

Table of Contents

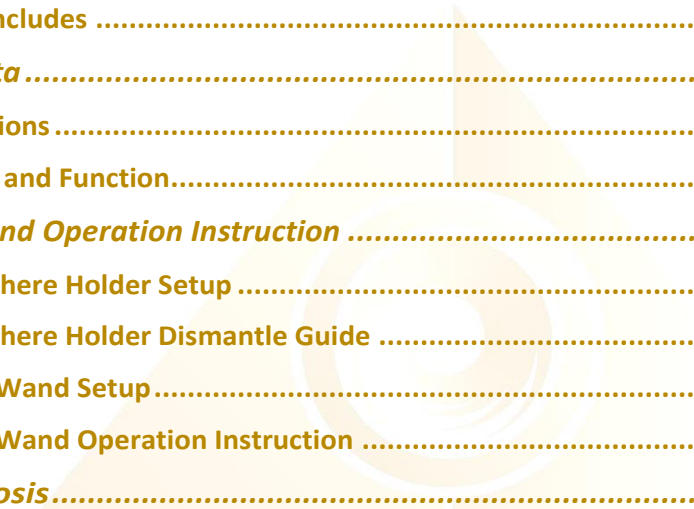

1	General Notes.....	2
1.1	Product Feature and Intend Use	2
1.2	Warning	2
1.3	Warranty	3
1.4	Precautions	3
1.5	Maintenance	4
1.6	Product Includes	5
2	Technical Data	6
2.1	Specifications	6
2.2	Overview and Function.....	7
3	Installation and Operation Instruction	9
3.1	Crystal Sphere Holder Setup	9
3.2	Crystal Sphere Holder Dismantle Guide	10
3.3	Mandala Wand Setup.....	11
3.4	Mandala Wand Operation Instruction	13
4	Fault Diagnosis.....	16

1 General Notes

1.1 Product Feature and Intend Use

Thank you for using Light Mandalas Product. Mandala Wand has been in good condition before leaving factory.

Mandala Wand is a portable device with **Non-Invasive Coherent Low Level Pulsed Lasers** (Red, Violet, and Near-Infrared), **Full Spectrum Light Emitter** (Photobiotic Lights), **Scalar Wave** Emitter (Shielding & Protection), **Healing Frequency** (1-40 MHz frequency), **Bifilar coil** as an energy converter, and **UVC LED** sterilization.

Mandala Wand has remarkable effects on human energy body such as:

- correct the energetic and informational fields from interfering or discordant energies and restores their balance
- mitigate the effects of stress, illness, pain and injury
- stimulate cell for healing and regeneration, improving appearance, performance, and overall well-being
- uplift human frequency to maintain health

1.2 Warning

Pay close attention to each safety REMARK and WARNING statement in the manual. This manual contains important laser system safety and operation information.

Read and understand all instructions prior to powering on the device the first time, to avoid laser eye injury and to avoid breaking the law. Lasers can be hazardous and have unique safety considerations. Permanent eye injury and blindness is possible if lasers are used incorrectly.

User should follow this manual to operate the product and make sure that everything is clear on all laser safety and product operation information before powering it on. Any wrong operating will make this unit out of warranty

This product is in perfect condition before delivery, please check all accessory when unpacking the unit.

Do not use in cases of Epilepsy. The flashing light may trigger a seizure pregnancy. Never look directly at light beam or its reflections. Follow standard laser eye protection.

1.3 Warranty

Each single product passes through a quality process. Nevertheless, if failures occur, contact the customer service at once. The product is under 12 months (1 year) warranty effectively according to the delivery date stated on the warranty card.

Warranty does not apply to damages, which result from misuse or neglect. The warranty also expires if you dismantle the product or failed to provide the warranty card. The manufacturer is not liable for consequential damage or in case of a non-intended use of the product.

If a failure occurs during the warranty period, the product will be replaced or repaired. The manufacturer reserves the right to exchange components of the product instead of repairing it. If the failure results from misuse or neglect, the user must pay for all cost that may occur. The freight cost and other costs will be depending on which Term & Conditions criteria the product falls into. **Import tax and Customs clearance duty cost is not included in the term of warranty.** In that case, you may ask for an estimated cost beforehand.

1.4 Precautions

1. Please check Mandala Wand for damage during shipping. Do not use the device if any damage is discovered. If the device is damaged, please contact the dealer or Light Mandalas immediately
2. Please Confirm the product you have received are in correct voltage and frequency for your country.
3. Do not disassemble or modify the device in any case.
4. Do not use the device outdoors when rains
5. Places the device in a well-ventilated area.
6. Keep dry and avoid excessive humidity, heat or dust.
7. Do not place the device closer than 2 feet from combustible object.
8. Do not look directly into the laser to avoid harm eyes.

9. Disconnect power when cleaning the tool or not in use.

1.5 Maintenance

1. Clean with soft cloth, humid tissue moistened with water or a water-based glass cleaner.
2. In order to guarantee the fixture works properly, it's essential to keep it clean.



1.6 Product Includes;

- | | |
|-------------------------------|----|
| 1. Mandala Wand GEN.2 | x1 |
| 2. Type-C Cable | x1 |
| 3. Charger 5V 2A | x1 |
| 4. Crystal Sphere Holder | x1 |
| 5. Crystal Sphere Set (30 mm) | |
| a. Alpha Sphere | x1 |
| b. Gamma Sphere | x1 |
| c. Theta Sphere | x1 |
| 6. Bag | x1 |



2 Technical Data


2.1 Specifications

Input:	Charger Adapter AC 100-240 V Type-C Cable Jack plug 3.5mm (audio output)
Main Memory:	SD Card
Main Frequency:	0-40MHz (square wave) Bifilar Coil 0-40MHz (sine wave)
Operation:	Bluetooth and Bluetooth Audio
Case Color:	Gold
Case Material:	Plastic ABS
Product size:	57 x 175 x 32 mm (W x H x D)
Gross Weight:	142g

2.2 Overview and Function

Front View and Back View

Front View

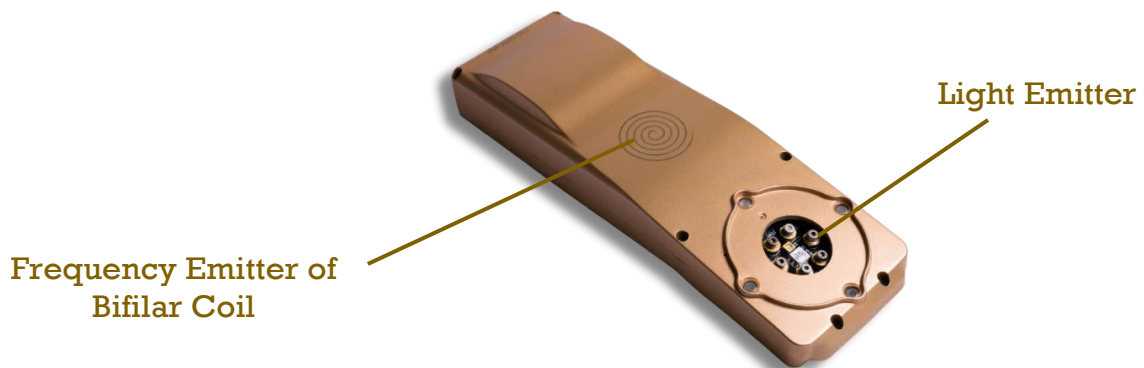

Side View and Bottom View

Side View



Power ON-OFF Button

Bottom View



Type-C Power Plug in 5V 2A

Audio Socket

3 Installation and Operation Instruction

3.1 Crystal Sphere Holder Setup

The holder is designed to fit perfectly on the Mandala Wand. Every holder has a mark on one of the four tips. The mark serves as a guide for the user to align and place the Mandala Wand's top magnet.



Red Circle - The Mark of the Crystal Sphere Holder

Align the mark to the top magnet of the Mandala Wand and place to stick it.



3.2 Crystal Sphere Holder Dismantle Guide

Place the holder downward and place your finger on the tip of the holder as shown in the pictures.



Gently pull out the holder. Make sure the holder is facing downward when dismantling. This is to avoid the crystal sphere drop out from the holder.



3.3 Mandala Wand Setup

1. Plug in Type-C cable to the power socket, the source can be USB Charger or Power Bank with 5V 2A or more.



2. Slide up the ON-OFF switch to ON




3. If a user wants to play music or setup pulsing laser, please insert the headset to the Mandala Wand Audio Socket



3.4 Mandala Wand Operation Instruction

When the device is switched ON, wait for a few seconds for the **Initializing** process.

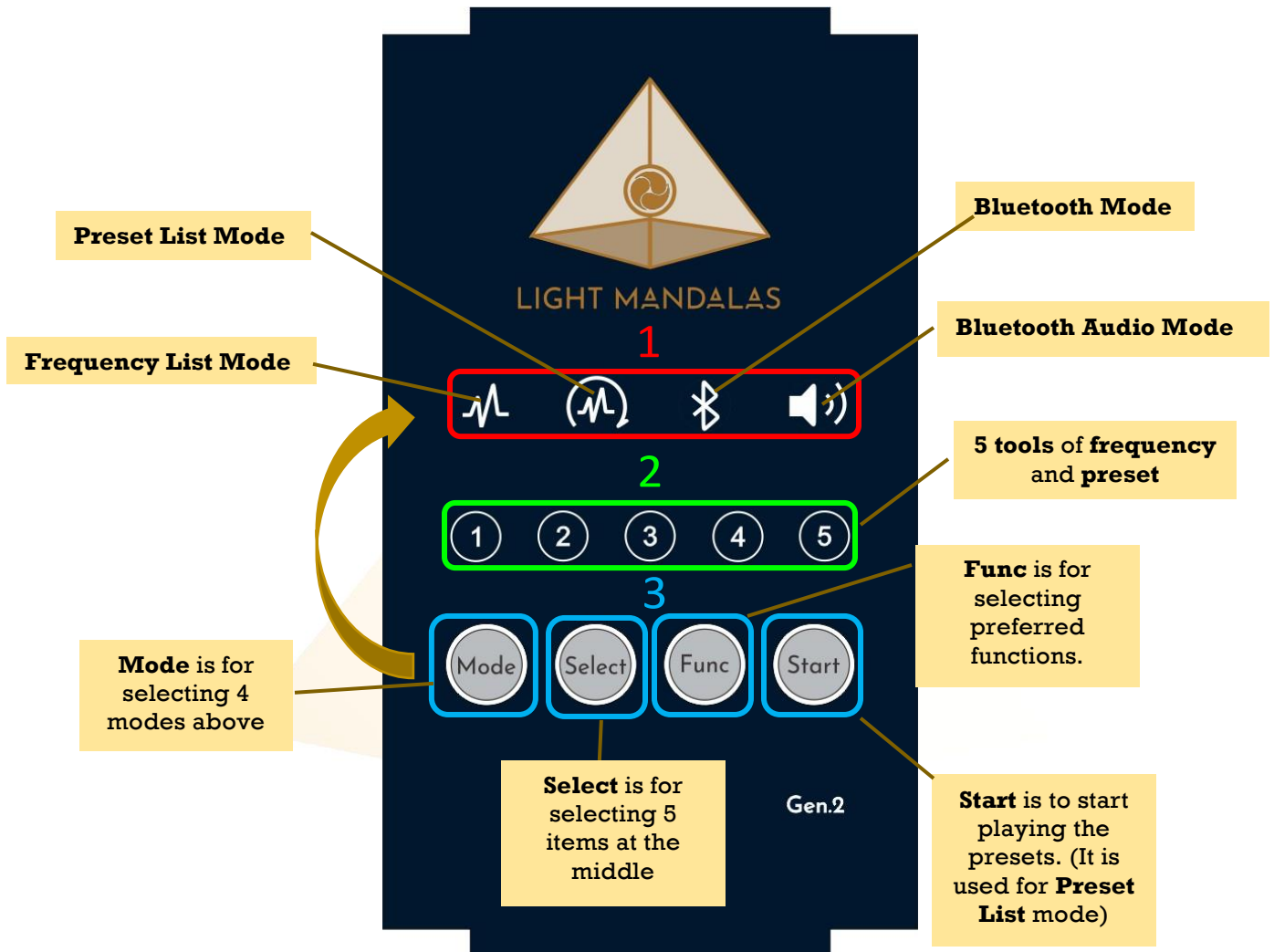

Mandala Wand Features:

Mandala Wand Operations:

1. The four circles at the top with a red light indicating the mode a user is in, including **Frequency List Mode**, **Preset List Mode**, **Bluetooth Mode**, and **Bluetooth Audio Mode**, which are all controlled by pressing the **Mode** button.
2. The five circles in the center, each with a green light, indicate **frequencies**

and **presets** (up to 5 items) that a user uploads to the Mandala Star App. They are activated by pressing the **Select** button.

3. The four buttons are used for controlling the usage.

3.1. **Mode** button: it is for selecting a preferred mode. There are four modes.

3.1.1. **Frequency List Mode**: In this mode, a user can select **specific frequency** to play. User can also upload specific frequencies to Mandala Wand from the **Mandala Star APP Frequency Library**, therefore, in the future user can play that specific frequency without using any phone or tablet.

3.1.2. **Preset List Mode**: Preset is a combination of several frequency lists. In this mode, a user can select **specific presets** to play. User can also manually create/customize their presets or download a **Basic Preset** shared by Light Mandalas in Mandala Star APP and upload to Mandala Wand. Therefore, in the future user can play that preset without using any phone or tablet.

**Please take note that you can only upload your preferred frequencies or presets as a set (maximum 5 items) at a time. When you upload a new set of frequencies/presets, your previous set of frequencies/presets will be automatically deleted.*

3.1.3. **Bluetooth Mode**: This mode is to allow a user to use Mandala Star App to control Mandala Wand. When this mode is active, the red light will flash.

Red light flashes
(Bluetooth mode)



3.1.4. **Bluetooth Audio Mode**: This is the combination of Bluetooth Mode and Audio Mode. In this mode, user can use Mandala Star App to control Mandala Wand and play music from the tablet. The lasers and scalar wave will pulse according to the music playing. When this mode is active, the red and blue light status will flash.

Red light flashes
(Bluetooth mode)



Blue light flashes
(Audio mode)



Red and Blue light flash
(Bluetooth Audio mode)



****Please take note that the Mandala Wand screen must switch to Bluetooth Audio mode and headset has to be plugged in in order to pair with tablet to play music.***

- 3.2. Select: it is for selecting a preferred frequency or preset.
- 3.3. Func: it is for selecting a preferred function.
- 3.4. Start: It is used for Preset List Mode. Tab on this button every time you start to play the presets.

4 Fault Diagnosis

Failures	Remedies
No power	<ul style="list-style-type: none">• Check the adapter input voltage and current. Does it match the required input voltage on the specification minimum 5V 2A?• Check On/Off switch. Make sure it is fully turn ON.• Check the Type-C cable. Make sure it is tightly attached.
1. Can't pair with tablet App 2. Can't upload frequency list or preset to Mandala Wand	<ul style="list-style-type: none">• Make sure the Mandala Wand screen is in Bluetooth Mode
3. Can't play music through tablet app	<ul style="list-style-type: none">• Make sure the Mandala Wand screen is in Bluetooth Audio Mode

Disclaimer

All content found on the Light-mandalas.com Website, Facebook Page, YouTube, including: text, images, video, audio, or other formats are intended for informational, educational and research purposes only and neither Light Mandalas Co. Ltd. nor the publishers of this content take responsibility for possible health consequences of any person or persons reading or following the information in this educational content and is not intended to be a substitute for professional medical advice, diagnosis, treatment or cure. Always seek the advice of your physician or other qualified health professionals with any questions you may have regarding a medical condition.

“No claims other than stress detection and stress reduction may be made”

Light Mandalas products are to be used as a Heart Rate Variability Biofeedback system and Light, Color and Sound therapy, designed for stress detection and stress reduction only. The devices do not evaluate any issue other than stress. Stress can come from many sources; this system uses light and color and sound therapies to reduce stress. Remember, there has never been any reported side effects using light, color and sound therapy anywhere in the world and throughout time.

Notes on the practices of low-level laser therapy (LLLT) and Biofeedback

Low-level laser therapy (LLLT) is considered investigational by the FDA, and no medical claims of cures are permitted in the US; at this time however, low-level lasers are defined within the scope of complimentary medical practices such as physiotherapy, aesthetic, Traditional Chinese Medicine and acupuncture, etc.

The FDA does not regulate Bio resonance tools but they do regulate biofeedback tools and some Bio resonance tools are registered under this category.

*It is a legal requirement at this time that the user does not make any medical claims of diagnosis or cure.

Precaution

Light Mandalas products are working on energetic level. The standard benefit is energy clearance. During the process, your body and mind can be uncomfortable. The symptoms associated with overly rapid clearance of toxins. If those detox symptoms are too strong, you may need to stop using the products. You have to allow the clearance or detox to proceed and you can start to try the products afterwards. However, if the issues still remain the same, it may need to be completely discontinued for your own safety

Shall there be any inquiry, please contact us for further support.

Light Mandalas Co., Ltd. (Thailand Headquarters):

- **Light Mandalas Website:** www.light-mandalas.com
- **Official:** info@light-mandalas.com
- **Instagram:** [@mandalawave9](https://www.instagram.com/mandalawave9)
- **YouTube Channel:** Light Mandalas Official
- **Technical Support Email:** cservice@light-mandalas.com

Local agents contacts:

- **Japan:** japan@light-mandalas.com
- **Taiwan:** lovelight777888@gmail.com
- **Indonesia:** lightseed2020bali@gmail.com
- **West Malaysia:** lightworker@healingtech.org
- **East Malaysia:** malaysia@light-mandalas.com
- **Hungary:** horesnyibela@light-mandalas.com
- **Switzerland, Germany, Austria, and the Netherlands:** info@light-mandalas.ch
- **France, Spain:** france-spain@light-mandalas.com
- **USA:** usa@light-mandalas.com
- **New Zealand, Australia:** tachyonhealing.nz@gmail.com